


|                     | 1 <sup>st</sup> APPOINTMENT   | 2 <sup>nd</sup> APPOINTMENT<br>TRY-IN  | 3 <sup>rd</sup> APPOINTMENT  |
|---------------------|---|--|--|
| DENTAL OFFICE STEPS | <ul style="list-style-type: none"> <li>• Take impression with wash inside the denture with border molding.</li> <li>• Bite registration.</li> <li>• Opposing impression model.</li> </ul>   | <ul style="list-style-type: none"> <li>• Verify vertical dimension.</li> <li>• Verify bite, and adjust occlusion if needed, then take a full bite registration.</li> <li>• Evaluate phonetics and esthetics.</li> </ul>                                | <ul style="list-style-type: none"> <li>• Seat the denture.</li> <li>• Verify bite, vertical dimension, esthetics, phonetics.</li> <li>• Adjust denture as needed.</li> </ul> |
| SEND TO PRODENCO    | <ul style="list-style-type: none"> <li>• Old denture with wash.</li> <li>• Bite registration.</li> <li>• Opposing model.</li> <li>• RX</li> </ul>   | <ul style="list-style-type: none"> <li>• Maxillary/mandibular models.</li> <li>• Try-in denture.</li> <li>• If new bite is required, take full arch bite registration, with denture try-in.</li> <li>• Additional instructions for changes.</li> </ul> |  |
| PRODENCO STEPS      | <ul style="list-style-type: none"> <li>• Create a new model.</li> <li>• Articulate to the bite.</li> <li>• Create a digital set up for try-in.</li> </ul>   | <ul style="list-style-type: none"> <li>• New articulation.</li> <li>• Digital set up.</li> <li>• Print Upper/Lower denture.</li> </ul>   |  |
| TIPS                | <ul style="list-style-type: none"> <li>• Verify tissue condition before impression (use tissue conditioner inside of old denture to heal the tissue and wait a couple days for impressions).</li> <li>• Verify vertical dimension when taking a new wash impression.</li> </ul> |  |  |

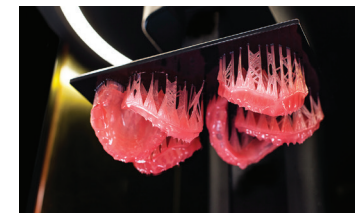

***Call us today  
with questions  
about this great  
new product!***

**PRODENCO.COM**  
611 Omaha St • Sioux City, Iowa  
**712.252.4034 • 800.831.0936**  
2218 West 2nd St • Grand Island, NE  
**308.282.0235 • 800.289.3525**  
407 N. 203 St • Elkhorn, NE  
**402.551.5954 • 800.747.5441**